

Chi Vitaliser


User Manual

Thank you for purchasing this product.


Please read this manual carefully to fully understand the proper usage of our product.

Product Details

Main Unit

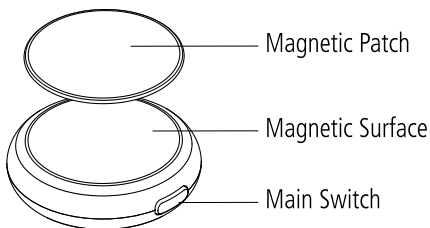

Main Unit Bottom View

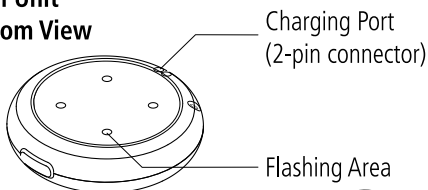


Charging Case

Charging Case Light Indicator



Connecting port between
the Case and the Main
Unit (2-pin connector)

Main Unit Light
Indicator

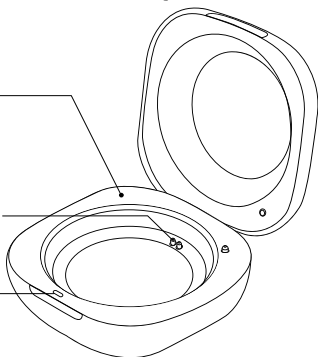

How to Use

- 1** Double tap to switch on/off

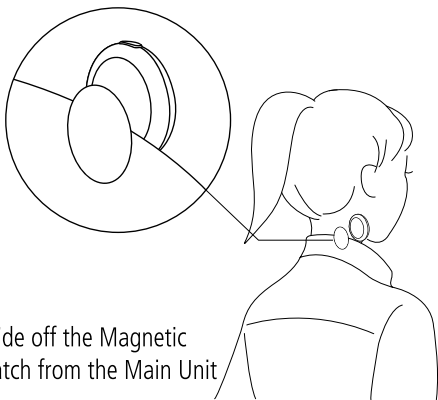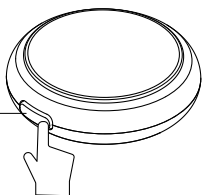


- 2** Slide off the Magnetic Patch from the Main Unit
- 3** Place the Main Unit over your body, flashing side towards your skin
- 4** Put the Magnetic Patch back onto the Main Unit, overlaying your clothes


Charging Instructions

A Charging Case
Light Indicator
(Red / Green)

B Main Unit
Light Indicator
(White)

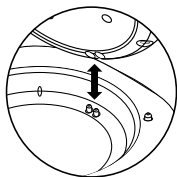



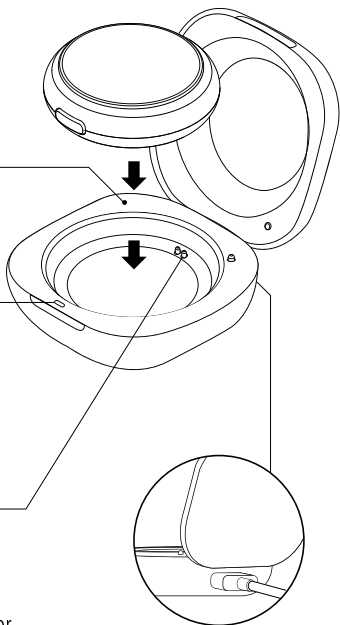

Align the 2-pin connector
at the bottom of the
Main Unit with the 2-pin
connector in the Case






Plug the Charging Cable
into the Type C port at
the back of the Case

A Charging Case Light Indicator

 Green Light remains lit	Fully charged
 Red Light remains lit	Charging
 Red Light blinks	Running out of battery*
 Light goes out	Running out of battery*

* It needs to be charged by using the Charging Cable. After connecting with the Charging Cable, the Red Light will become lit.

B Main Unit Light Indicator

 White Light remains lit	Fully charged
 White Light blinks slowly	Charging
 Light goes out	Running out of battery#

The Charging Case is low in battery to charge the Main Unit. It needs to be charged by using the Charging Cable.


After connecting with the Charging Cable for 5-8 seconds, the White Light will blink slowly to indicate charging is in progress. If there is no light on, please check if the 2-pin connector at the bottom of the Main Unit and the 2-pin connector in the Case are aligned.

Anti-allergy Sticker (Suitable for sensitive skin)


How to Use




Open the Sticker package

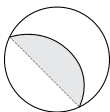

Peel off the self-adhesive back paper

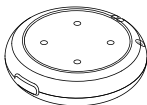

Clean the base of the Main Unit and fix the Sticker onto the base

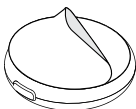

Change the Sticker weekly to keep it clean

Caution

1. Do not put this product in hot, humid or dusty environment, or place it under direct sunlight.
2. Do not subject the product to strong impact to avoid malfunction.
3. Do not disassemble or modify the product to avoid danger.
4. Keep the product out of reach of children.
5. Stop using the product in the event of damage or malfunction.

Product Maintenance

1. Keep the product clean and free from acidic, oily, sticky or other contaminants.
2. If the product is contaminated, wipe it clean with a soft, dry cloth.
3. When in use, please avoid excessive collision with hard objects.
4. Avoid liquids, such as water, from flowing into the Main Unit.

Packing List

Main Unit	x 1	Anti-allergy Sticker	x 4
Charging Case	x 1	User Manual	x 1
USB Type C Charging Cable	x 1	Acupoint Manual	x 1
Magnetic Patch	x 2		

Product Specification

Product Name	: Chi Vitaliser
Model No.	: EF19205
Dimension	: 14.6 x ø 46.2 mm (Main Unit) 28.3 x 57.8 x 57.8mm (Charging Case)
Weight	: 30g (Main Unit) 38g (Charging Case)
Materials	: ABS + PC
Battery	: Li-Ion 3.7V 420 mAh (Main Unit) Li-Ion 3.7V 560 mAh (Charging Case)
Battery Life	: 6 - 8 hours (Main Unit) 2 - 3 hours (Charging Case)
Charging Time	: 4 hours (Main Unit) 4 hours (Charging Case)
Operating Voltage	: 4.2V (Li-Ion battery)
Rated Input	: DC 5V 2A
Power Source	: 4 x 0.02W / LED
Standard	: GB 4706.15-2003 / GB 7247.1-2012
Product Name	: Anti-allergy Sticker
Dimension	: ø 34mm



Customer Service



If you have any further queries, please do not hesitate to contact us



We will do our best to solve the issue

Hotline

(852) 2677 7898

Manufacturer: New Age Technology (Asia) Limited

Address : Unit B6, 35/F, C&D Wings, TML Tower,
3 Hoi Shing Rd, Tsuen Wan, Hong Kong

Office hours: Monday – Friday 09:00 – 18:00

(Lunch break 13:00-14:00)

Closed on Saturdays, Sundays, Public Holidays