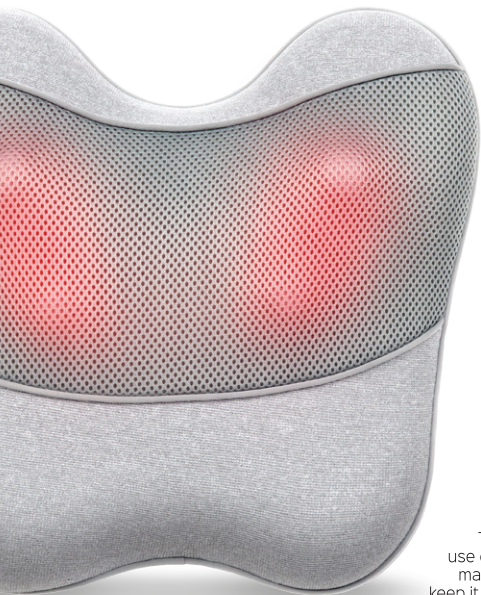


Energy Waist Massager

User Guide



To ensure that you get the best use of this product, please read this manual thoroughly before use and keep it handy when using the product.

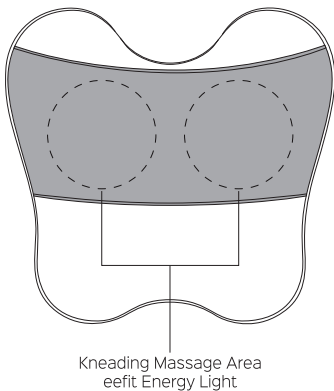
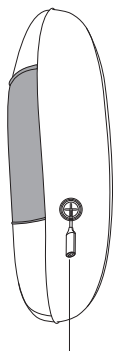
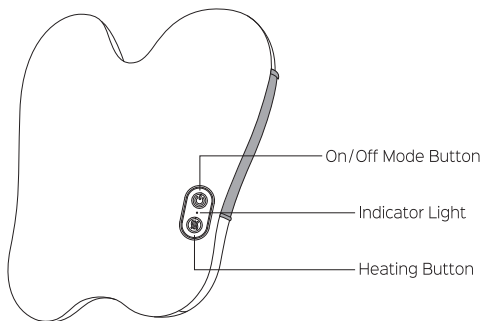
Disclaimer:
The contents of this manual and the specifications of this product are subject to change without notice. The company reserves the right to make changes without notice in the specifications and materials contained herein and shall not be responsible for any damages (including consequential) caused by reliance on the materials presented, including but not limited to typographical and other errors relating to the publication.

Please use the product in accordance with the below instructions.

1. Please use a Type-C power supply compatible with this product. If the case or power cord is damaged, do not use the product.
2. To avoid damage, do not alter or replace the cord or circuitry of this product, or pull the power cord or plugs too hard.
3. Please keep the product clean and dry, do not use the product or touch the power supply if your hands are wet.
4. Please unplug the power supply when not in use.
5. In the event of a power outage, stop using and disconnect the power supply immediately to avoid possible product damage when power is restored.
6. Do not insert any objects into the crevices of the massager.
7. Do not use accessories and attachments other than those specified.
8. It is recommended for use indoors and in safe environments.
9. Do not squeeze the product too hard to avoid injury.
10. It is recommended to wait one hour after meal before use.
11. If you experience any discomfort or abnormality during use, please stop using immediately and consult a doctor or professional.
12. If you experience any abnormal reactions during use, stop using immediately and consult the supplier.




Features

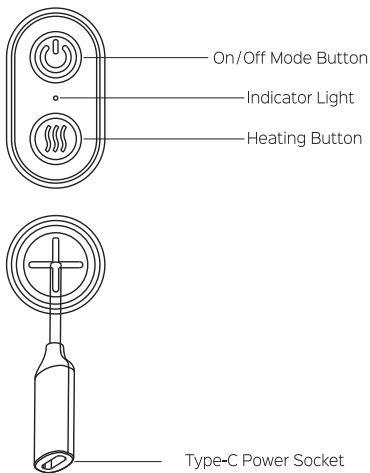
Product structure



Operation Instructions

When using for the first time or if battery power is lost from not using for a long time, please charge the product according to the instruction manual to ensure it can work normally.

1. On-Off Button: Long press  to turn the machine on and off , when the green light is on, it is in Mode 1.
2. MODE Button: The default mode is Mode 1, short press  to switch to Mode 2.
3. Heating Button: The default mode has no heating. Short press  to turn on/off heating. The red light at the massage head indicates that the heating is on.
4. Indicator Light: When in use, the green light will be on. When the power is low, the red light will flash until the power runs out. When charging: the red indicator light will turn on. Fully charged: the green indicator light will turn on. (Please don't use when charging)
5. Type-C Charging Port: Connect the Type-C charging cable to the charging port to start charging.



Caution

1. This product is designed as an electronic healthcare massage product for use at home and is not suitable for diagnosing or treating a disease.
2. Do not use while sleeping.
3. Never use it with the following medical electronic instruments:
 - Implantable medical electronic instruments such as cardiac pacemakers.
 - Life-sustaining medical electronic devices such as artificial heart-lung machines.
4. The following users should consult a physician before using this product. Otherwise, accidents or physical discomfort may occur:
 - Patients with abnormal heart and lumbar nerves and who cannot state their wishes.
 - Patients with severe hyperosteogeny.
 - People with abnormal spine or back curvature caused by accidents or diseases.
 - Pregnant women, infants and children, patients with skin sensory disorders or skin abnormalities.
5. When not in use, please put the product and the power supply in the packaging and store it in a dry and ventilated place to avoid moisture.
6. If the product is not functioning properly, please contact us. Please read the warranty terms and conditions for more information.

Product Specification

Product Name:	Energy Waist Massager
Model:	P035
Size:	320 x 130 x 305mm
Weight:	1.85kg
Battery Capacity:	2000mAh
Operation Time:	10mins
Fixed input voltage:	5V
Fixed input power:	10W
Heating Temperature :	40-50°C



Troubleshooting

	Problem	Possible causes	What to do
1	Cannot be switched on	i. Insufficient power ii. Did not press and hold the On/Off Button	i. Use when fully charged ii. Press and hold the On/Off button
2	Unable to charge or no indicator light will not turn on	The Type-C charging cable is not connected to the charging port	Connect the Type-C charging cable to the charging port
3	The massage suddenly stops working	i. Insufficient power ii. The automatic 10-minutes timer has ended	i. Use when fully charged ii. If you would like to continue massaging press the On/Off button again

Please scan the QR code
to activate your product warranty



Warranty Terms and Conditions

- The manufacturer offers an exchange within 7 days from the date of purchase if there is any product quality problem.
- The manufacturer provides free spare parts replacement and repair service for any defects found in the product hardware within 1 year from the date of purchase.
- The manufacturer does not provide any warranty for the damage that may result from the loss of spare parts, intentional or unintentional misuse of the product, self-disassembly of the product or modification of alteration of the product by unauthorized personnel.
- The manufacturer does not provide any warranty for damage caused by accidents (including fire, water damage, contact with any liquid, etc.).
- Please retain the receipt and present it to the staff as proof if repair service is required.
- The manufacturer reserves the right to amend the above terms and conditions at any time without prior notice to the customer.

Hotline: (852) 2677 7898